

# STRADIS<sub>MED</sub>



Blood Transfusion

## Docon Seal M

Portable Tube Sealer



STRADIS-MED GmbH  
Boschring 12  
D-63329 Egelsbach  
Amtsgericht Offenbach HRB 42891  
Phone: +49 (0) 6103 73333-0  
Fax: +49 (0) 6103 73333-10  
E-Mail: [info@stradis-med.de](mailto:info@stradis-med.de)  
Web: [www.stradis-med.de](http://www.stradis-med.de)

## Docon Seal M

### HIGH PERFORMANCE SEALING WITH THE Docon Seal M

The Docon Seal M is an automated sealing unit to seal PVC tubes, particularly the blood removal tube directly attached to the donor. This ensures increased level of safety when sealing the tube. Its handy, robust and light construction combined with a long battery charge life provide the best conditions for everyday use in production and further processing of blood products.

Different tube diameters and wall thicknesses can be sealed.

The sealing time is automatically adjusted to the tubes specified. An LED indicator on the sealing handle shows that the sealing handle is in process.

### KEY FEATURES


- Light, robust and handy mobile sealing unit
- Ergonomically shaped sealing handle including spray protection (also suitable for small hands)
- Can be operated using minimal manual force
- Easy to clean
- Latest lithium-ion battery technology for up to 2000 seals per battery charge
- Reliable system for sealing directly at the donor
- Automatically regulated sealing time depending on tube thickness

### TECHNICAL DATA

Dimensions of handle (LxWxH) 200x25x70mm • Dimensions of device (LxWxH) 160x90x55mm

Weight of handle 350g • Weight of device approx. 950 g • Tube specification 4.0 - 6.0 mm

Sealing time 1.5 s depending on tube diameter

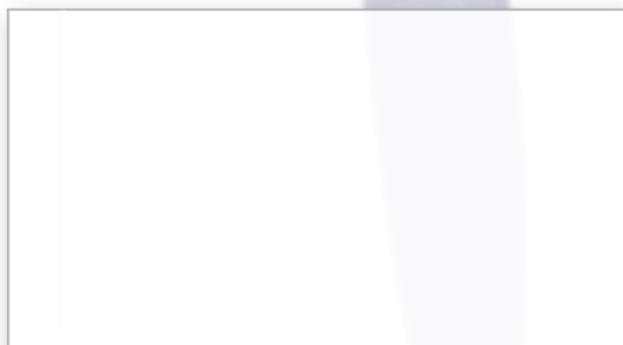

**Möller Medical GmbH**  
Wasserkuppenstraße 29-31  
36043 Fulda, Germany  
Phone +49 66194195-0  
Fax +49 66194195-850  
info@moeller-medical.com  
www.moeller-medical.com

Distributed by **STRADIS<sub>MED</sub>**

**STRADIS-MED GmbH**  
Boschring 12  
D-63329 Egelsbach  
Amtsgericht Offenbach HRB 42891  
Phone: +49 (0) 6103 73333-0  
Fax: +49 (0) 6103 73333-10  
E-Mail: info@stradis-med.de  
Web: www.stradis-med.de