

# STRADIS<sub>MED</sub>

## RAPIDWELD™ STW

- Waferless Sterile Welding
- Touch Screen Interface
- Infrared Heating
- Unique Patented Clamps



### The RapidWeld™ STW

- yields accurate, reliable, sterile welds without the requirement of a heated consumable copper wafer
- can improve operator efficiency and maximize the shelf life of valuable blood products
- sterile tube welder is an approved sterile connection device with the ability to produce non-contaminating welds
- is easy to use, cost effective and an eco-friendly device designed to produce sterile welds using an innovative approach to sterile tube welding



**GENESIS BPS** EC REP

Genesis BPS International sàrl  
CP 2361 Nyon 2, Switzerland  
Phone: + 41.22.363.18.16  
Fax: + 41.22.363.18.17  
international@genesibps.com

STRADIS-MED GmbH  
Woogstrasse 48  
D-63329 Egelsbach  
Amtsgericht Offenbach HRB 42891  
Office: +49 6103 70 64 22  
Fax: +49 6103 70 64 229  
info@stradis-med.de  
www.stradis-med.de

## The RapidWeld™ STW Sterile Connection Advanced Technology Modern Design

- Waferless Sterile Welding
- Touch Screen Interface
- Infrared Heating
- Unique Patented Clamps

### The RapidWeld™ STW

- yields accurate, reliable, sterile welds without the requirement of a heated consumable copper wafer
- can improve operator efficiency and maximize the shelf life of valuable blood products
- sterile tube welder is an approved sterile connection device with the ability to produce non-contaminating welds
- is easy to use, cost effective and an eco-friendly device designed to produce sterile welds using an innovative approach to sterile tube welding

### Features and Benefits

#### Advanced Blade Module

The stainless steel blade permits continuous and uninterrupted welding with no jamming.

- Increase operator efficiency and speed
- Notifies user when blade replacement is required
- Blade module kit: 1000 welds

#### User Interface with LCD Touch Screen Display

LCD color touch screen (4" x 3") display provides operator with step-by-step instructions for performing welds.

- Easy to use menu driven interface with audio & visual alarms
- Languages: English, French, German, Italian, Portuguese & Spanish

#### Infrared Heating Element

- Designed for continuous operation
- Maintains a high temperature to assure weld sterility

#### Lid Lock Mechanism

- Protects operator during welding process
- Avoids unintentional interruption of the weld

#### Patented Tube Holder Clamps

Unique clamp design strips fluid away from the weld site and keeps tubing securely in place.





This process can help eliminate blood cell necrosis in blood products.

- Parallel tube loading for quick and simple use
- Clearly marked indicators for operator guidance

## CREATE STRONG RELIABLE WELDS IN 14-15 SECONDS

Fast and easy step-by-step welding process

CONNECT ANY COMBINATION OF WET AND DRY PVC TUBING

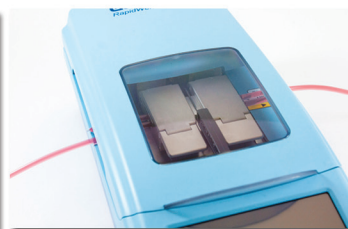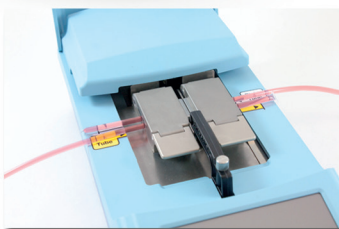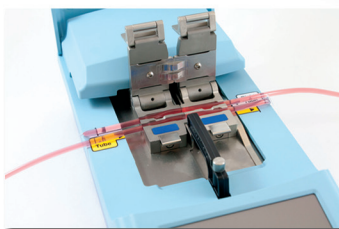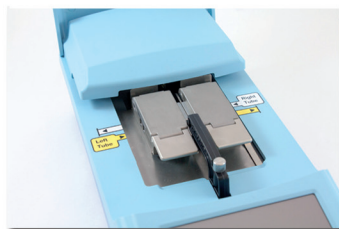

Applications include, but are not limited to:

- Pathogen reduction processing - Component aliquoting - Quality control sampling - Bacterial testing - Leukoreduction - Collection set modification
- Cell washing and freezing - Apheresis set modification - Transferring research samples - Component pooling - Cellular Therapy - Other PVC applications

Article	Model	Article Code	
Genesis Rapid Weld	GRW-430	480-RW-2607-00	- Warranty: 1 year limited warranty
Blade Modul Kit	1000 Welds	480-RW-2478-03	- FDA and CE proofed